

# TEG<sup>®</sup> Functional Fibrinogen

Whole-blood fibrinogen monitoring assay



TEG<sup>®</sup> functional fibrinogen test helps guide the most informed treatment decisions possible

INTEGRATED DEVICES

**TEG<sup>®</sup> Functional Fibrinogen assay determines the contribution of functional fibrinogen to clot strength enabling targeted and informed treatment decisions.**


### **TEG<sup>®</sup> 5000 - Function Fibrinogen Solutions:**

Fast and easy-to-use

Clinically relevant information to help assess the hemostatic condition:

- Thrombotic or bleeding risk
- Ratio of platelets to fibrinogen
- Cause of bleeding

Optimal measurement and sensitivity of hyper-fibrinolyses

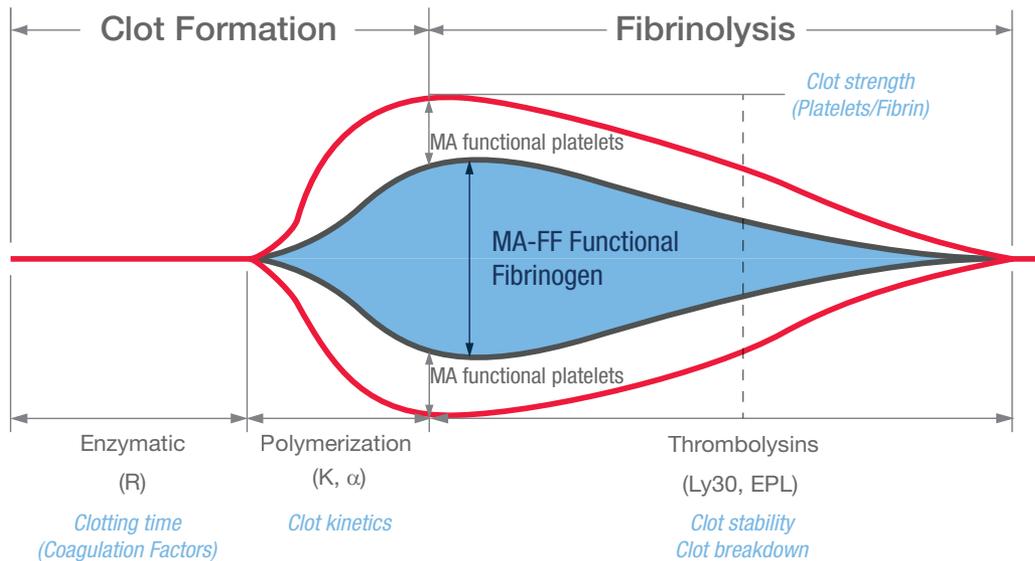

## Tests available

- Functional Fibrinogen (non-citrated)
- Citrated Functional Fibrinogen
- Functional Fibrinogen Heparinase (non-citrated)
- Citrated Functional Fibrinogen Heparinase

## Ordering information

Description:	List Number:
Functional Fibrinogen Test	07-034

## Related products

Haemonetics has a complete offering of patient blood management solutions. To order or for technical support, please contact your regional office or local representative.

## Supporting information

The functional fibrinogen test measures the fibrinogen that actually contributes to the clot strength (Functional Fibrinogen). The platelet contribution is blocked by a potent inhibitor and the resulting MA-FF reflects just the fibrin contribution to the clot strength. The strength of the fibrin clot is influenced mainly by fibrinogen- and F XIII levels of the blood sample and by fibrin polymerization disorders.



CE 0123

© 2013 Haemonetics Corporation. Haemonetics and TEG are trademarks or registered trademarks of Haemonetics Corporation in the USA, other countries, or both. 07.2013, Germany. COL-COPY-000686(AA)

Haemonetics S.A. ■ Signy Centre ■ Rue des Fléchères ■ P.O. Box 262 ■ 1274 Signy-Centre ■ Switzerland ■ www.haemonetics.com

**Austria**  
0800.29.2777

**Belgium and Luxembourg**  
FR: 0800.754.80  
NL: 0800.754.82

**Czech Republic**  
800.143.243

**Denmark**  
8088.7112

**France**  
0800.90.11.58

**Germany**  
0800.180.8890

**Italy**  
800.870.200

**The Netherlands**  
0800.0222.707

**Norway**  
800.18.453

**Sweden**  
020.797.150

**Switzerland**  
0800.898.898

**United Kingdom**  
0808.234.4817